INVENTUM

SINCE 1908

BabyCare

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PAGE 8
PAGE 8
PAGE 8
PAGE 9
PAGE 9
PAGE 9
PAGE 14

ENGLISH



Legislation requires that all electrical and electronic equipment must be collected for reuse and recycling. Electrical and electronic equipment marked with the symbol indicating separate collection of such equipment must be returned to a municipal waste collection point.

DESCRIPTION OF THE APPLIANCE

1. SAFETY INSTRUCTIONS
 2. BEFORE FIRST USE
 3. WEIGH
 4. BATTERIES
 5. CLEANING & MAINTENANCE
- GUARANTEE & AFTER SALES SERVICE

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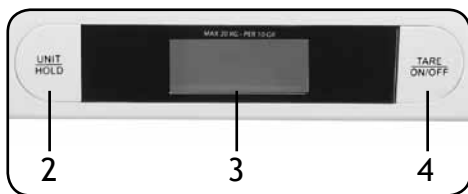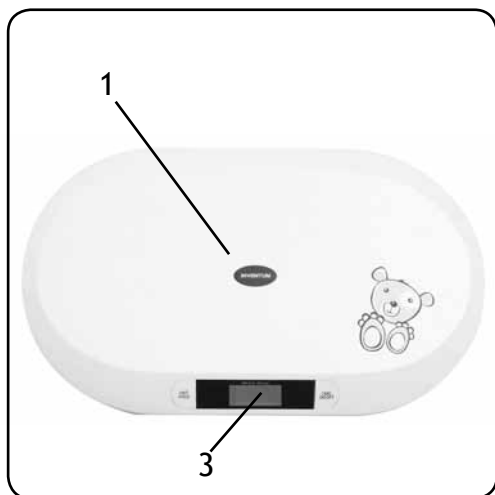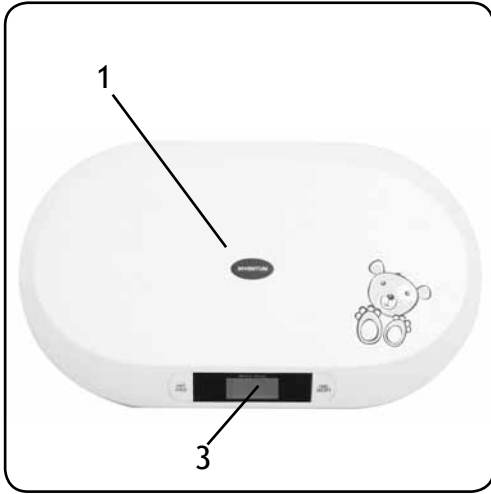
1. Weighing surface
2. UNIT/HOLD key
3. Display
4. ON/OFF/TARE key

1 | SAFETY INSTRUCTIONS

- **Please read these instructions before operating the appliance and retain for future use.**
- The appliance is intended solely for domestic use within the home.
- Always use batteries which corresponds with the voltage printed on the bottom of the appliance.
- Always use the appliance on a level and solid surface; this is important for ensuring that the scale weighs correctly.
- For your own safety, only use the accessories and spare parts from the manufacturer which are suitable for the appliance.
- Protect your scale from impact with hard objects, moisture, dust, chemicals, toiletries, liquid cosmetics, great temperature fluctuation and closeness to sources of heat (open fires, radiators).
- Strong electromagnetic fields (e.g. cell phones) may impair the accuracy of the scale.
- Always return the appliance after a malfunction, or if it has been damaged in any manner to the Customer Service Department for examination, repair or adjustment as special purpose tools are required.
- An electrical appliance is not a toy and should therefore always be placed outside of the reach of children. Children are not always able to understand potential risks. Teach children a responsible handling of electrical appliances.
- Children should be supervised to ensure that they do not play with the appliance.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- If you decide never to use the appliance again, we advise you to dispose of it in special containers or at a

municipal waste collection point.

- Never use the appliance if it is not working properly, or if it has been damaged. Any intervention other than cleaning and normal maintenance by the customer must be carried out by an approved service centre.
- **Never immerse any part of the main body in water or any other liquid to protect against electrical hazards.**
- If the appliance is used for other purposes than intended, or is it not handled in accordance with the instruction manual, the full responsibility for any consequences will rest with the user. Any damages to the product or other things are not covered by the warranty.

2 | BEFORE FIRST USE

Before you use the appliance for the first time, please do as follows: Carefully unpack the appliance and remove all the packaging material. Keep the material (plastic bags and cardboard) out of reach of children. Check the appliance after unpacking for any damage, possible from transportation. Place the enclosed batteries - 2 x 1,5Volt AAA - in the batterie compartment on the bottom of the baby scale. Wipe the outside of the baby scale only with a damp cloth.

3 | WEIGH

The baby scale has a capacity of minimum 10 grams to a maximum of 20 kilograms (44 lb/3.2 st). Weight is displayed in 10 g graduations.
Put the appliance on a plane and stable surface.

Switch the scale on by pushing the **ON/OFF/TARE** key. Wait about 2 seconds until «**0.00 kg**» appears on the display. The baby scale is now ready for weighing. Now place the baby on the weighing surface. Be extremely careful not to let the baby fall and take care of the baby's head.
The display will show the weight of the baby.

This baby scale has a **HOLD** function. This function fixes the weight of the baby as soon as it is still enough to record an accurate measured value. Press **HOLD** and the display of the baby scale will blink while displaying the weight for 20 seconds as soon as the baby lies still enough. The scale then automatically switches off. If the baby moves a lot, it is impossible to determine its weight. You should therefore try to pacify the baby if it is kicking in order to do the measurement promptly. Do not hold the baby, however, as otherwise this will affect the weight reading.

If you do not want to place the baby directly on the weighing surface, but prefer to use a towel, for example, you have to place the towel on the scale **BEFORE** switching the scale on. If the towel and any other items which you want to deduct again during weighing, you can also use the **TARA** function. Press the **TARA** button once the items to be deducted have been weighed. Place the baby on the scale and continue as described above.

To select a different measurement unit press the **UNIT** key to change it accordingly to: kg = kilograms, lb = pounds and st = stone.
Remove the overweight when «**EEEE**» appears in the display. **NOTE:** The total weight placed on the scale can not be more than 20 kg.

4 | BATTERIES

If the baby scale is left unused for longer than approx. 20 seconds, it automatically switches itself off. You can also turn the scale off manually by touching the key **ON/OFF/TARE** for about 2 seconds.

When the  symbol appears in the display, the batteries need to be replaced. Open the battery compartment at the bottom of the scale. If necessary please use a small screwdriver to remove the batteries. Make sure the batteries are placed correct (= and - polarity) and close the lid.

NOTE: Spent batteries and rechargeable batteries do not constitute normal household waste! They are considered to be toxic waste and, as such, should be disposed of in special containers, toxic waste collection points or brought to electrical goods dealers.

Remark: Batteries containing pollutant substances are marked as follows: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury.

5 | CLEANING & MAINTENANCE

You can clean the baby scale with a damp cloth and a little washing up liquid, if required.

Never immerse the scale in water or any other liquids. Handle the scale with care, protect the scale from impact with hard object and do not drop it.

Do not place any objects on the scale when not in use.

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GUARANTEE & AFTER SALES SERVICE

- This product is guaranteed for a period of 2 years from the date of purchase against mechanical and electrical defect.
- Guarantee can only be given if proof of purchase (e.g. till receipt) can be shown. Therefore keep it in a safe place. The guarantee period starts at date of purchase.
- Guarantee is not valid in case of:
 - normal wear;
 - damage through misuse or when the instructions were not followed correctly;
 - lack of maintenance;
 - not following the operation- and safety instructions;
 - incompetent montage or repairs made by yourself or others;
 - using accessories and spare parts which are not suitable for the appliance; for assessment by our technical department.
- The guarantee is not valid for:
 - normal wear;
 - batteries;
 - (interchangeable) signal- and controllamps for assessment by our technical department.

Return to your supplier or send the appliance, including a copy proof of purchase and description of the complaint, stamped to:

Inventum Group BV
Generatorstraat 17
3903 LH VEENENDAAL
The Netherlands
Website: www.inventum.eu

IES Import Electronic Service GmbH
Lembekstrasse 36
22529 Hamburg
GERMANY
Website: www.ihrservice.de

- Pack it carefully (preferably in the original carton). Ensure unit is clean.
- Enclose your name and address and quote model number on all correspondence.

When the guarantee is no longer valid or when the period of guarantee has passed, we will charge any costs that flow out of an assessment and/or repairing.

On request we will supply a statement before repairing. Without requesting a statement, we assume that you, by sending the appliance, automatically order to repair the appliance. A guarantee period of two months is valid for repairs done by our technical department and the used spare parts, counting from the repair date. When a malfunction occurs within this period which directly flows out of the repairs and the used spare parts, the consumer has the right of repairs without any costs.

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GL-International GmbH